



oil baths

WWW.FALCINSTRUMENTS.IT

oil baths





I N D E X

04

oil baths

Oil baths
with round tank

05

Oil baths
with rectangular tank

oil baths

OIL BATHS WITH ROUND TANK



BO



BO-D



Excellent heat insulation guaranteed by the double wall. The tanks are made of 18/8 stainless steel sheet, no welds and seals, pressed in one piece



Heating is ensured by electrical heating elements fixed under the bottom of the tanks, thus the tanks are completely free and easy to clean



BO: Temperature is controlled by a thermostat with sensor.

BO-D: The temperature The temperature is set and checked by a digital thermostat with PID controller and PT100 probe.



The instrument is fitted with an adjustable safety thermostat

Technical data		BO-4	BO-12	BO-D24	BO-D30
Capacity	lt	6	13	6	13
Heating power	W	1000	1400	1000	1400
Precision	°C	5	5	1	1
Drain		yes	yes	yes	yes
Temperature	°C	+50 to +300	+50 to +300	+5 amb. to +300	+5 amb. to +300
External dimensions LxWxH	mm	305x260	390x340	305x290	390x340
Internal dimensions LxWxH	mm	220x150	250x200	220x150	280x200
Weight	kg	5,2	8,2	5,8	9,4
Code		614.1033.03	614.1033.12	614.1035.04	614.1035.12

V/HZ 220-240/50-60

SAFETY CLASS3,1

PROTECTION CLASS IP 42

ACCESSORIES



18/8 Stainless steel lid with concentric rings

BO-4/BO-D24

Code 610.1010.24

BO-12/BO-D30

Code 610.1010.32



18/8 Stainless steel test tube rack with 3 shelves with central handle

	BO-4/BO-D24	BO-12/BO-D30
	36 holes Ø21	60 holes Ø21
Code	212.9260.03	212.9260.04

-  The external case and tank are made of 18/8 stainless steel sheet, no welds and seals, pressed in one piece
-  Manual reset safety thermostat with alarm light, red coloured
-  PID temperature controller with PT100 probe



WB-OD


Technical data		WB-OD5	WB-OD15	WB-OD24
Tank capacity	lt	5	12	18,5
Adjustable temperature	°C	+5 amb. to +300	+5 amb. to +300	+5 amb. to +300
Precision in tank H ₂ O	°C	±1.0	±1.0	±1.0
Drain tap		yes	yes	yes
Heating power	W	700	1400	2100
Internal dimensions LxWxH	mm	100x250x150	255x290x150	255x440x150
External dimensions LxWxH	mm	280x360x240	445x385x240	445x560x240
Weight	kg	5,2	6,2	9
Code		614.1032.05	614.1032.15	614.1032.24

V/HZ 220-240/50-60

SAFETY CLASS 1

PROTECTION CLASS 42

ACCESSORIES

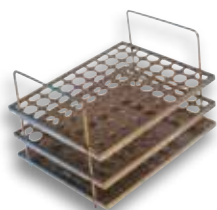

18/8 Stainless steel test tube rack with 3 shelves

	2 per tank	1 per tank	2 per tank
	12 holes Ø21	48 holes Ø21	48 holes Ø21
Code	212.1250.02	212.1250.05	212.1250.05
	20 holes Ø17	60 holes Ø17	60 holes Ø17
Code	212.1250.03	212.1250.06	212.1250.06



18/8 Stainless steel lid with concentric rings. Max per hole Ø85 mm

Fori n.	2	4	6
Code	610.1007.35	610.1007.34	10.1007.36



18/8 stainless steel lid

Code	616.1025.05	616.1025.12	616.1025.14
-------------	--------------------	--------------------	--------------------

Notes

[illegible]

Notes

[illegible]



FALC INSTRUMENTS s.r.l.

Via G. M. Compagnoni, 2
24047 Treviglio (BG) - Italy
+39 0363 304660 | falc@falcinstruments.it

WWW.FALCINSTRUMENTS.IT